



PHILIPS

VitalEye


MR Systems

Keeps a caring eye on your patient

With VitalEye, the technologist no longer needs to set up an old-fashioned respiratory belt, but receives a continuous and robust respiratory signal without any interaction. This revolution in patient sensing helps you keep a caring eye on your patient and provides superior image quality¹, even for challenging patients, thanks to touchless patient sensing. The quality of the physiology signal is significantly better than a belt-based approach based upon faster, more accurate, detection of breathing.

VitalEye

System	Ingenia Elition X
Coil compatibility	dS Torso, dS HeadSpine
Patient orientation	Head-first / Feet-first
Body locations	Up to 50
Frame rate	20 frames /s
Temporal resolution	50 ms
Infra red compatibility	Yes
Independent of exam room light	Yes
Motion detection resolution	Sub-mm
Accessories needed for respiratory triggering	No
Calibration	Automatic

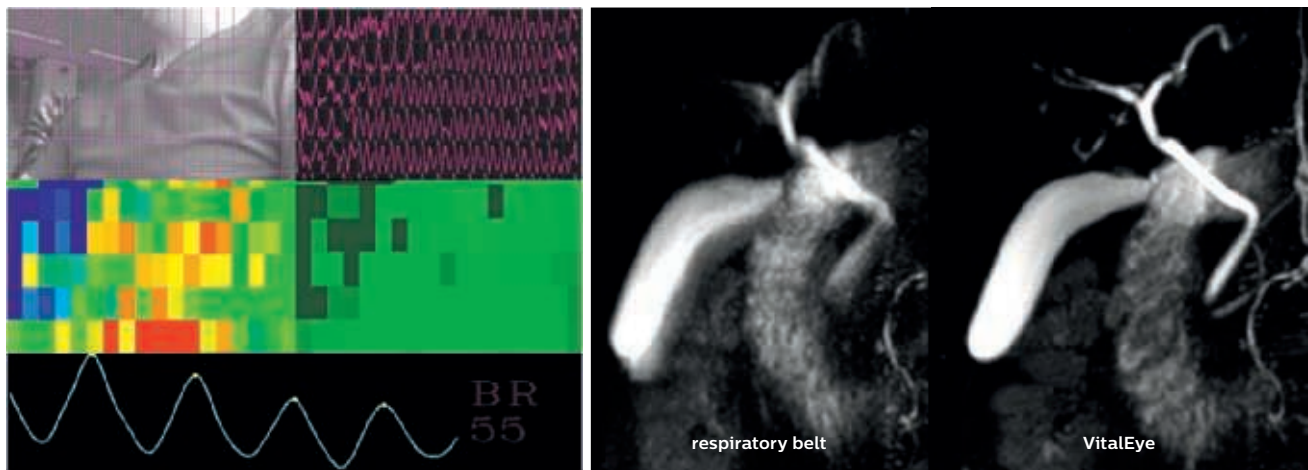

VitalEye uses adaptive intelligence to detect even the tiniest microscopic signs of breathing. Advanced algorithms analyze up to 50 body locations parallel and in real time. This produces a respiratory trace that is more robust compared to traditional respiratory triggering methods.

1 Compared to respiratory belt